


actiTENS

A wearable and app-controlled TENS for chronic pain management.

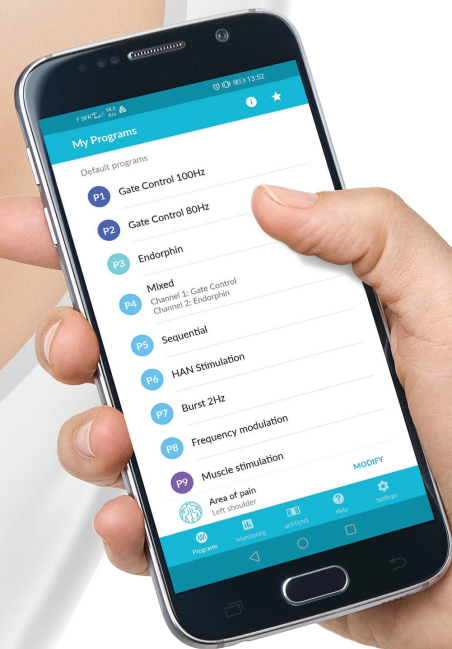


SUBLIMED

MADE IN FRANCE



From CHRONIC PAIN to ACTIVE LIFE...

20% of the world's population suffers from chronic pain.

Sublimed has worked with patients and healthcare providers to offer a new generation of drug-free devices to relieve patients from chronic pain.

actiTENS, more than 15 000 users.

actiTENS is a W-TENS - Wearable Transcutaneous Electrical Nerve Stimulator - providing a clinically proven therapy* for chronic pain management.

* Maheu E. et al., (2022). Wearable transcutaneous electrical nerve stimulation (actiTENS®) is effective and safe for the treatment of knee osteoarthritis pain. A randomized controlled trial versus weak opioids. *Therapeutic Advances in Musculoskeletal Diseases*, volume 14



Unique low-back electrode

RELIEVE YOUR PAIN **UNNOTICED**


actiTENS discreetly accompanies patients in their daily life, allowing them to resume an active life.

actiTENS can stimulate every part of the body. Its easy-to-wear, compact, thin, and flexible form adapts itself perfectly to each patient's morphology.

★★★★★ **Bérangère. D, 03/03/2021**

“ Using Actitens has changed my daily life. Significant reduction of chronic pain. Easy to use, it is worn without constraints. No matter what I do, I can do my sessions. ”

...with actiTENS

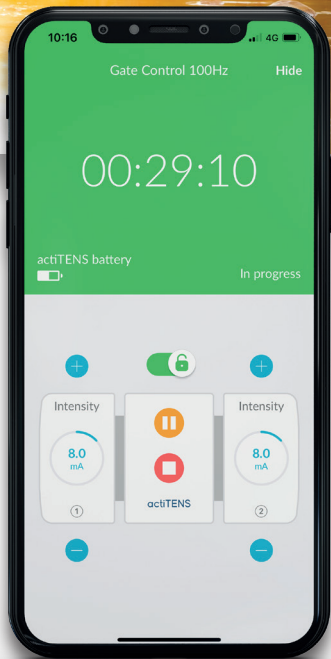

START A STIMULATION PROGRAM IN 15 SECONDS

Bluetooth-enabled **actiTENS** intuitive app is available on



Google Play 4,5 ★



App Store 4,5 ★




Mohamed E,
28/04/2021

“ Perfect! Without a doubt a device that will relieve you of your aches and pains. ”

KEEP TRACK OF YOUR PROGRESS

Stimulation history is automatically recorded and is available under monitoring tab.

Pre- and post-session pain level can be monitored to track progress and engage patients in their recovery.

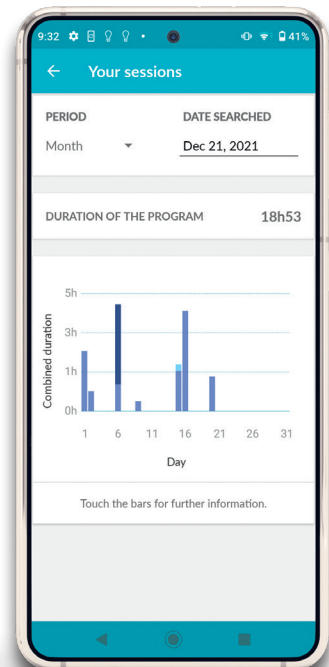

PROGRAMS

- P1 Gate Control 100 Hz
- P2 Gate Control 80 Hz
- P3 Endorphin
- P4 Mixed
- P5 Sequential
- P6 HAN Stimulation
- P7 2 Hz Burst
- P8 Frequency Modulation
- P9 Muscle stimulation
- P10 Massage
- P11 Friction
- P12 Sensitive areas
- P13 Nausea
- P14 Tibial nerve stimulation

TECHNICAL DATA

Channels: 2 independent channels

Battery: Li-Ion fully charged within 2 hours

Battery autonomy: 20 hours* (Gate Control 100 Hz, 2 channels, 15 mA)

Output voltage: 60 mA \pm 10% max. (1000 Ω) / 60 V max.

actiTENS dimensions:
108 mm x 53.5 mm x 17 mm

actiTENS weight: 64 g

Warranty: 2 years warranty on **actiTENS** and **actiCHARGE** only.

actiTENS IP rating: IP22

Pulse width: 50-400 μ s \pm 5 μ s

Frequency: 1-120 Hz \pm 10%

Treatment time: 10-720 minutes, adjustable via the mobile app.

Cables: 4 lengths available (14 cm, 40 cm, 70 cm and 100 cm)

Electrodes: wide range for all indications
*patient's skin impedance 500 Ω

- ① AC charger / ② Armband / ③ Self-adhesive strip / ④ actiTENS
⑤ Electrodes / ⑥ Cables / ⑦ actiCHARGE / ⑧ Transport case



actiTENS is a medical device subject to contraindications.
Please carefully read the **actiTENS** manual before use.


SUBLIMED

137, rue Mayoussard
38430 Moirans • France
+33 (0)4 76 37 17 58
contact@subli-med.com

CE
2797

FDA CLEARED
510(k) K202159

 **French**
Healthcare
Association member

www.subli-med.com